

Iron Stain (Liver) - Tissue Control Slides

Description: Human.
Natural disease state (where applicable).
Paraffin embedded.

Availability: TCS0013-5 5 unstained slides.
TCS0013-25 25 unstained slides.

Storage: 2-25° C.

Suggested Stain Kit: IRN-1

Procedure: See below.

Iron Stain Kit


Description: The Iron Stain Kit is intended for use in the detection of ferric iron in tissues, blood smears, or bone marrow smears. Ferric iron is normally found in small amounts in bone marrow and the spleen. Abnormally large deposits may be seen in hemochromatosis and hemosiderosis. This product is based on the Prussian Blue reaction in which ionic iron reacts with acid ferrocyanide producing a blue color.

Tissue Sections


Iron:	Bright Blue
Nuclei:	Red
Background:	Pink


Uses/Limitations: Not to be taken internally.
For In-Vitro Diagnostic use only.
Histological applications.
Do not use if reagents become cloudy.
Do not use past expiration date.
Use caution when handling reagents.
Non-Sterile.

Control Tissue: Liver
Spleen
Bone Marrow

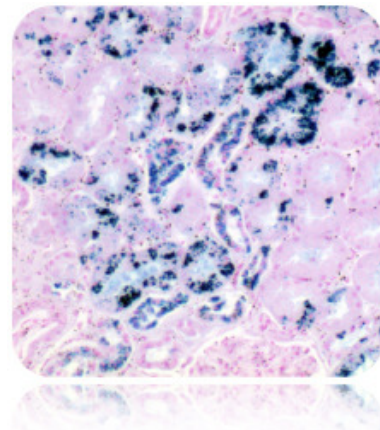


Ordering information regarding individual components on back page!

Storage: 2° C  25° C

 ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.

 EmergoEurope (31)(0) 70 345-8570
Molsnstraat 15
2513 BH Hague, The Netherlands

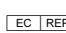Kit Contents:

<u>Item #</u>	<u>Kit Contents</u>	<u>Volume</u>	<u>Storage</u>
PFB500	Potassium Ferrocyanide Solution	500 ml	18-25°C
HQB500	Hydrochloric Acid Solution (2%)	500 ml	18-25°C
NFS125	Nuclear Fast Red Solution	125 ml	18-25°C

Precautions: Avoid contact with skin and eyes.
 Harmful if swallowed.
 Follow all Federal, State, and local regulations regarding disposal.

Procedure:

Note: Use acid-washed or bleach-washed glassware.
 Rinse all glassware with distilled water prior to use.
 Do not use metal forceps to transfer slide during staining procedure.


1. Deparaffinize sections if necessary and hydrate to distilled water.
2. Mix equal volumes of Potassium Ferrocyanide Solution and Hydrochloric Acid Solution to make a working Iron Stain Solution. Use once and discard.
3. Incubate slide in working Iron Stain Solution for 3-5 minutes.
4. Rinse thoroughly in distilled water.
5. Stain slide in Nuclear Fast Red Solution for 5 minutes.
6. Rinse in 4 changes of distilled water.
7. Dehydrate in 3 changes of absolute alcohol for 2 minutes each.
8. Clear, and mount in synthetic resin.


References:

1. Sheenan, D.C., Hrapchak, B.B. Theory and Practice of Histotechnology, 2nd Edition. Battelle Press, Columbus, OH. Page 217. 1980
2. Carson, F.L., Histotechnology; A Self-Instructional Text, ASCP Press, Chicago, IL. Pages 214-215. 1990

Bulk Reagent Ordering Information and Current Pricing at www.scytek.com


Description:	Catalog #	Volume
Potassium Ferrocyanide Solution	PFB500	500 ml
	PFB999	1000 ml
Hydrochloric Acid Solution (2%)	HQB500	500 ml
	HQB999	1000 ml
Nuclear Fast Red Solution (Enhanced Stability)	NFS125	125 ml
	NFS500	500 ml
	NFS999	1000 ml

 Storage: 2° C  25° C



ScyTek Laboratories, Inc.
 205 South 600 West
 Logan, UT 84321
 U.S.A.

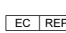

 EmergoEurope (31)(0) 70 345-8570
 Molsnstraat 15
 2513 BH Hague, The Netherlands