


Calcium Stain Kit (Modified Von Kossa)

Description: The Calcium Stain Kit (Modified Von Kossa) is intended for use in the histological visualization of calcium deposits in paraffin sections.

Calcium in mass deposits:	Black
Calcium in dispersed deposits:	Gray
Nuclei:	Red
Cytoplasm:	Light Pink

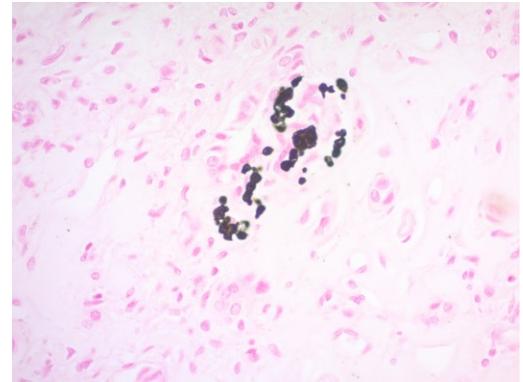

Uses/Limitations: For In-Vitro Diagnostic use only.
 Histological applications.
 Do not use past expiration date.
 Use caution when handling these reagents.

Control Tissue: Any paraffin embedded tissue that contains calcium deposits.

Availability/Contents:

<u>Item #</u>	<u>Kit Contents</u>	<u>Volume</u>	<u>Storage</u>
SNV030	Silver Nitrate Solution (5%)	30 ml	2-8°C
STB030	Sodium Thiosulfate Solution (5%)	30 ml	18-25°C
NFS030	Nuclear Fast Red Solution	30 ml	18-25°C

Precautions: Avoid contact with skin and eyes.
 Harmful if swallowed.
 Follow all Federal, State, and local regulations regarding disposal.
 Use in chemical fume hood whenever possible.

Procedure (Standard):

1. Deparaffinize sections if necessary and hydrate to distilled water.
2. Incubate slide in 6-8 drops of Silver Nitrate Solution (5%) for 30-60 minutes while exposing to either ultraviolet light or incandescent light at 75 watts or greater. For best results, keep light source within 2 feet (61cm) of slide during Silver Nitrate staining procedure.
3. Rinse in 3 changes of distilled water.
4. Incubate slide in 5-8 drops of Sodium Thiosulfate Solution (5%) for 2 minutes.
5. Rinse for 2 minutes in running tap water followed by 2 changes of distilled water.
6. Stain tissue section with 6-8 drops of Nuclear Fast Red Solution for 5 minutes.
7. Rinse for 2 minutes in running tap water followed by 2 changes of distilled water.
8. Dehydrate very quickly in 3 changes of Absolute Alcohol.
9. Clear, and mount in synthetic resin.

Storage: 2° C  25° C

**Mixed Storage Conditions.
 Separate Contents.**

Doc: IFU-TemplateMixedStoragerev2



ScyTek Laboratories, Inc.
 205 South 600 West
 Logan, UT 84321
 U.S.A.




EMERGO
 Emργο Europe
 Prinsessegracht 20
 2514 AP The Hague, The Netherlands

References:

1. Sheenan, D.C., Hrapchak, B.B. Theory and Practice of Histotechnology, 2nd Edition. Battelle Press, Columbus, OH.
2. Clark, G., et al. Staining Procedures, 4th Edition, Williams & Wilkins Press, Baltimore, MD.
3. Symonds, D.A., Use of the Von Kossa stain in identifying occult calcifications in breast biopsies. American Journal of Clinical Pathology, 1990, July; 94(1) pages 44-48.

Storage: 2° C  25° C**Mixed Storage Conditions.
Separate Contents.**

Doc: IFU-TemplateMixedStorageev2

 ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.CE 

EC REP

Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands