

Picro-Sirius Red Stain Kit (For Collagen)

Description: The Picro-Sirius Red Stain Kit (For Collagen) is intended for use in the histological visualization of muscle and collagen fibers in tissue sections. The PSR stain may be viewed using standard light microscopy or polarized light resulting in birefringence of collagen fibers.

Light Microscopy


Collagen:	Red
Muscle Fibers:	Yellow
Cytoplasm:	Yellow

Polarized Light Microscopy

Collagen:	Yellow-Orange and Green Birefringence
-----------	---------------------------------------

Uses/Limitations:

- Not to be taken internally.
- For In-Vitro Diagnostic use only.
- Histological applications.
- Do not use if reagents become cloudy.
- This protocol has not been tested on and may not be suitable for frozen sections.
- Do not use past expiration date.
- Use caution when handling reagents.
- Non-Sterile.

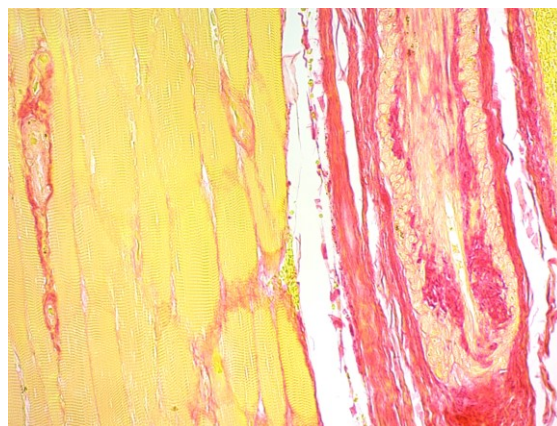

Control Tissue:

- Lung
- Uterus
- Muscle
- Kidney

Availability/Contents:

<u>Item #</u>	<u>Kit Contents</u>	<u>Volume</u>	<u>Storage</u>
SRS030	Picro-Sirius Red Solution	30 ml	18-25°C
AAD030	Acetic Acid Solution (0.5%)	30 ml x 2	18-25°C


Ordering information regarding individual components on back page!

Precautions:

- Avoid contact with skin and eyes.
- Harmful if swallowed.
- Follow all Federal, State, and local regulations regarding disposal.
- Use in chemical fume hood whenever possible.

Procedure:

Storage: 18° C  25° C

 ScyTek Laboratories, Inc.
 205 South 600 West
 Logan, UT 84321
 U.S.A.

CE 

EC REP
 Emurgo Europe
 Prinsessegracht 20
 2514 AP The Hague, The Netherlands

1. Deparaffinize sections if necessary and hydrate to distilled water.
2. Apply adequate Picro-Sirius Red Solution to completely cover tissue section and incubate for 60 minutes.
3. Rinse slide quickly in two changes of Acetic Acid Solution (0.5%).
4. Rinse slide using absolute alcohol.
5. Dehydrate in 2 changes of absolute alcohol, clear, and mount in synthetic resin.

References:

1. Puchtler H., Waldrop F.S., Valentine L.S. Polarization microscopic studies of connective tissue stained with picro-sirius red FBA. Beitr Path. 1973; 150, pages 174-187.
2. Junqueira L.C.U., Bignolas G., Brentani R.R. Picrosirius staining plus polarization microscopy, a specific method for collagen detection in tissue sections. Histochemistry J. 1979, 11, pages 447-455.
3. Whittaker P. Polarized light microscopy in biomedical research. Microscopy and Analysis 1995; 44, pages 15-17.

Bulk Reagent Ordering Information and Current Pricing at www.scytek.com

Description:	Catalog #	Volume
Picro-Sirius Red Solution	SRS250	250 ml
	SRS500	500 ml
	SRS999	1000 ml
Acetic Acid Solution	AAD250	250 ml
	AAD500	500 ml
	AAD999	1000 ml

Storage: 18° C  25° C



ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.



EC REP

Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands