

# SensiTek HRP Anti-Rabbit (DAB) Ready-To-Use (70 slide)


<b>Species of Origin:</b>	Goat
<b>Antigen Specificity:</b>	Anti-Rabbit IgG (H+L)
<b>Preabsorbed Against:</b>	Human
<b>Enzyme Conjugate:</b>	Peroxidase
<b>Chromogen Substrate:</b>	Diaminobenzidine


## Procedure:

1. Deparaffinize and rehydrate tissue section.
2. To reduce nonspecific background staining due to endogenous peroxidase, incubate slide in hydrogen peroxide for 10-15 minutes.
3. Wash 2 times in buffer.
4. If required, incubate tissue in digestive enzyme.
5. Wash 4 times in buffer.
6. Apply Super Block (blue cap), and incubate for 5-10 minutes at room temperature to block nonspecific background staining. *Note: Do not exceed 10 minutes or there may be a reduction in desired stain.*
7. Wash 1 time in buffer.
8. Apply primary antibody and incubate according to manufacturer's protocol.
9. Wash 4 times in buffer.
10. Apply Biotinylated Link Antibody (yellow cap), and incubate for 15-20 minutes at room temperature.
11. Wash 4 times in buffer.
12. Apply Streptavidin/HRP Label (red cap), and incubate for 20 minutes at room temperature.
13. Rinse 4 times in buffer.

**WARNING:** DAB is a suspected carcinogen. Handle with care and dispose of according to all regulations.

14. Add 8 drops of DAB chromogen to one 5ml vial of DAB Substrate. Mix well and apply to tissue for 5 minutes.
15. Rinse 1 time in DI Water.

Storage: 2° C  8° C



ScyTek Laboratories, Inc.  
205 South 600 West  
Logan, UT 84321  
U.S.A.

CE 

EC REP

Emergo Europe  
Prinsessegracht 20  
2514 AP The Hague, The Netherlands

P.O. Box 3286 - Logan, Utah 84323, U.S.A. - Tel. (800) 729-8350 - Fax (435) 755-0015 - [www.scytek.com](http://www.scytek.com)

16. Apply DAB Chromogen/Substrate mixture and incubate for a second 5-minute period.
17. Rinse 3 times DI water.
18. Counterstain and coverslip.

### **- Troubleshooting Guide -**

#### **Overstaining:**


1. Concentration of the primary antibody was too high or the incubation time was too long.
2. Temperature during incubation was too high.
3. Incubation time with link antibody or streptavidin/enzyme label was too long.


#### **Nonspecific Background Staining:**

1. Rinsing between steps was inadequate.
2. Tissue was allowed to dry with reagents on.
3. Folds in tissue trapped reagents.
4. Tissue contains endogenous peroxidase.
5. Tissue contains endogenous biotin.
6. Antigen migrated in tissue.
7. Excessive tissue adhesive on slides.
8. Inadequate blocking with protein block.

#### **Weak Staining:**


1. Primary antibody concentration was too low or incubation time was too short.
2. Reagents are past their expiration date.
3. Inadequate removal of wash water between steps, resulting in dilution of reagents.
4. Room temperature was excessively cool.
5. The primary antibody does not recognize an antigen that survives fixation and embedding in high enough amounts.
6. Excessive incubation with protein block (Super Block).

Storage: 2° C  8° C



ScyTek Laboratories, Inc.  
205 South 600 West  
Logan, UT 84321  
U.S.A.

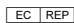


  
Emergo Europe  
Prinsessegracht 20  
2514 AP The Hague, The Netherlands

P.O. Box 3286 - Logan, Utah 84323, U.S.A. - Tel. (800) 729-8350 - Fax (435) 755-0015 - [www.scytek.com](http://www.scytek.com)

**No Staining:**

1. Steps were inadvertently left out.
2. There is no antigen in the tissue.
3. The primary antibody is not of rabbit origin.
4. Chromogenic substrate has been replaced with another that is not intended for use with peroxidase.

Storage: 2° C  8° C



ScyTek Laboratories, Inc.  
205 South 600 West  
Logan, UT 84321  
U.S.A.

CE 



Emergo Europe  
Prinsessegracht 20  
2514 AP The Hague, The Netherlands