



Instructions For Use

PSR-IFU

205 South 600 West Logan, Utah 84323, U.S.A. – Tel. (800) 729-8350 – Tel. (435) 755-9848 – Fax (435) 755-0015 – www.scytek.com Rev. 5, 7/25/2022

Picro-Sirius Red Stain Kit (For Collagen)

Description and Principle

The Picro-Sirius Red Stain Kit (For Collagen) is intended for use in the histological visualization of muscle and collagen fibers in tissue sections. The PSR stain may be viewed using standard light microscopy or polarized light resulting in birefringence of collagen fibers.

Collagen and reticular fibers are selectively stained by Sirius Red in a strongly acidic solution saturated with picric acid. Sirius Red binds to collagen in a highly ordered parallel manner enhancing the natural birefringent properties of collagen.

Expected Results

Light Microscopy

Collagen:	Red
Muscle Fibers:	Yellow
Cytoplasm:	Yellow

Polarized Light Microscopy

Collagen: Yellow-Orange and Green Birefringence

Kit Contents

1. Picro-Sirius Red Solution
2. Acetic Acid Solution (0.5%)

Storage

- 18-25°C
- 18-25°C

Suggested Controls (not provided)

Lung, Uterus, Muscle, Kidney.

Uses/Limitations

For In-Vitro Diagnostic use only.

Do not use if reagents become cloudy or precipitate

Do not use past expiration date.

Use caution when handling reagents.

Non-Sterile

Intended for FFPE sections cut at 5-10µm.

This procedure has not been optimized for frozen sections.

Frozen sections may require protocol modification.

Storage


Store kit and all components at room temperature (18-25°C).

Safety and Precautions

Please see current Safety Data Sheets (SDS) for this product and components GHS classification, pictograms, and full hazard/precautionary statements.

Procedure:

1. Deparaffinize sections if necessary and hydrate to distilled water.
2. Apply adequate Picro-Sirius Red Solution to completely cover tissue section and incubate for 60 minutes.
3. Rinse slide quickly in two changes of Acetic Acid Solution (0.5%).
4. Rinse slide using absolute alcohol.
5. Dehydrate in 2 changes of absolute alcohol, clear, and mount in synthetic resin.

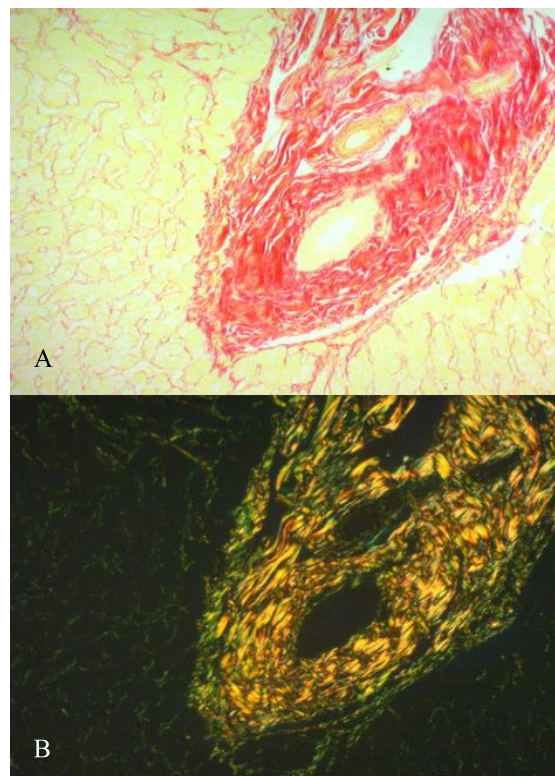

Collagen in human liver demonstrated with Picro-Sirius Red and viewed under A) Brightfield, and B) Polarized light.

References

1. Ferdousy, Raihana Nasrin, Hiroya Kadokawa, and Patrick Lonergan. "Specific locations and amounts of denatured collagen and collagen-specific chaperone HSP47 in the oviducts and uteri of old cows as compared with those of heifers." *Reproduction, Fertility and Development* (2022). <https://doi.org/10.1071/RD21130>.
2. Song, J.W., Ahn, J.W., Lee, M.W. et al. Targeted theranostic photoactivation on atherosclerosis. *J Nanobiotechnol* 19, 338 (2021). <https://doi.org/10.1186/s12951-021-01084-z>
3. Jiheun Ryu, Ungyo Kang, Joon Woo Song, Junyoung Kim, Jin Won Kim, Hongki Yoo, and Bomi Gweon, "Multimodal microscopy for the simultaneous visualization of five different imaging modalities using a single light source," *Biomed. Opt. Express* 12, 5452-5469 (2021). <https://doi.org/10.1364/BOE.430677>.
4. Shang, Q., Chu, Y., Li, Y. et al. Adipose-derived mesenchymal stromal cells promote corneal wound healing by accelerating the clearance of neutrophils in cornea. *Cell Death Dis* 11, 707 (2020). <https://doi.org/10.1038/s41419-020-02914-y>
5. Mengting Zhu, Yunpeng Chu, Qianwen Shang, Zhiyuan Zheng, Yanan Li, Lijuan Cao, Yongjing Chen, Jianchang Cao, Oscar K. Lee, Ying Wang, Gerry Melino, Guozhong Lv, Changshun Shao, Yufang Shi, Mesenchymal stromal cells pretreated with pro-inflammatory cytokines promote skin wound healing through VEGFC-mediated angiogenesis, *Stem Cells Translational Medicine*, Volume 9, Issue 10, October 2020, Pages 1218–1232, <https://doi.org/10.1002/sctm.19-0241>
6. Ogino, R.; Hayashida, K.; Yamakawa, S.; Morita, E. Adipose-Derived Stem Cells Promote Intussusceptive Lymphangiogenesis by Restricting Dermal Fibrosis in Irradiated Tissue of Mice. *Int. J. Mol. Sci.* 2020, 21, 3885. <https://doi.org/10.3390/ijms21113885>
7. Okano S, Yasui A, Kanno SI, Satoh K, Igarashi M, Nakajima O. Karyopherin Alpha 2-Expressing Pancreatic Duct Glands and Intra-Islet Ducts in Aged Diabetic C414A-Mutant-CRY1 Transgenic Mice. *Journal of diabetes research.* 2019;2019.
8. Moriya C, Imai K, Taniguchi H. PRDM 14 is overexpressed in chronic pancreatitis prior to pancreatic cancer. *FEBS open bio.* 2018 Oct;8(10):1733-41.

9. Narikawa M, Umemura M, Tanaka R, Fujita T, Yokoyama U, Ishigami T, Kimura K, Tamura K, Ishikawa Y. Acute Hyperthermia Inhibits TGF- β 1-induced Cardiac Fibroblast Activation via Suppression of Akt Signaling. *Scientific reports*. 2018 Apr 19;8(1):6277.
10. R. A. Helmer, R. Martínez-Zaguilán, J. S. Dertien, C. Fulford, O. Foreman, V. Peiris, and B. S. Chilton, "Helicase-Like Transcription Factor (Hltf) Regulates G2/M Transition, Wt1/Gata4/Hif-1 α Cardiac Transcription Networks, and Collagen Biogenesis," *PLoS ONE*, vol. 8, no. 11, p. e80461, Nov. 2013.
11. Puchtler H., Waldrop F.S., Valentine L.S. Polarization microscopic studies of connective tissue stained with picro-sirius red FBA. *Beitr Path*. 1973; 150, pages 174-187.
12. Junqueira L.C.U., Bignolas G., Brentani R.R. Picrosirius staining plus polarization microscopy, a specific method for collagen detection in tissue sections. *Histochemistry J*. 1979, 11, pages 447-455.
13. Whittaker P. Polarized light microscopy in biomedical research. *Microscopy and Analysis* 1995; 44, pages 15-17.



ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
435-755-9848
U.S.A.



Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands