


CD38 (ADP Ribosyl Cyclase I); Clone CD38/4328 (Concentrate)

Availability/Contents:	Item #	Volume
	RA0731-C.1	0.1 ml
	RA0731-C.5	0.5 ml
	RA0731-C1	1 ml

Description:

Species:	Mouse
Immunogen:	Recombinant full-length human CD38 protein
Clone:	CD38/4328
Isotype:	IgG1 / Kappa
Entrez Gene ID:	952
Hu Chromosome Loc.:	4p15
Synonyms:	ADP-ribosyl cyclase/cyclic ADP-ribose hydrolase 1, 2'-phospho-ADP-ribosyl cyclase, 2'-phospho-ADP-ribosyl cyclase/2'-phospho-cyclic-ADP-ribose transferase, 2'-phospho-cyclic-ADP-ribose transferase, ADP-ribosyl cyclase 1, Cyclic ADP-ribose hydrolase 1, T10, Acute Lymphoblastic Leukemia Cells Antigen; ADP Ribosyl Cyclase 1; cADP-ribose Hydrolase 1; CD38H; NAD(+) Nucleosidase; NIM-R5 Antigen; p45; T10
Mol. Weight of Antigen:	~45kDa (Glycoprotein); 35kDa (protein core)
Format:	200ug/ml of antibody purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide.
Specificity:	CD38 is expressed on CD34+ cells.
Background:	CD38 is a type II transmembrane glycoprotein that is present on early B- and T-cell lineages and activated B- and T-cells but is absent from most mature resting peripheral lymphocytes. CD38 is also found on thymocytes, pre-B cells, germinal center B-cells, mitogen-activated T-cells, monocytes and Ig-secreting plasma cells. The CD34+CD38- population of hematopoietic stems cells define the most pluripotent cells (e.g. blast colony forming cells).
Species Reactivity:	Human
Positive Control:	Raji cells. Human tonsil, Spleen or Skin.
Cellular Localization:	Membrane
Titer/ Working Dilution:	Immunohistochemistry (Frozen and Formalin-fixed): 1-2µg/ml Flow Cytometry: 1-2µg/million cells Immunofluorescence: 1-3µg/ml Western Blotting: 2-4µg/ml
Microbiological State:	This product is not sterile.s

Storage: 2° C  8° C



ScyTek Laboratories, Inc.
 205 South 600 West
 Logan, UT 84321
 U.S.A.

Instructions For Use RA0731-C-IFU-RUO

Rev. Date: Aug 20, 2024


Revision: 1

Page 2 of 2

P.O. Box 3286 - Logan, Utah 84323, U.S.A. - Tel. (800) 729-8350 – Tel. (435) 755-9848 - Fax (435) 755-0015 - www.scytek.com

Uses/Limitations:

- Not to be taken internally.
- For Research Use Only.
- This product is intended for qualitative immunohistochemistry with normal and neoplastic formalin-fixed, paraffin-embedded tissue sections, to be viewed by light microscopy.
- Do not use if reagent becomes cloudy.
- Do not use past expiration date.
- Non-Sterile.

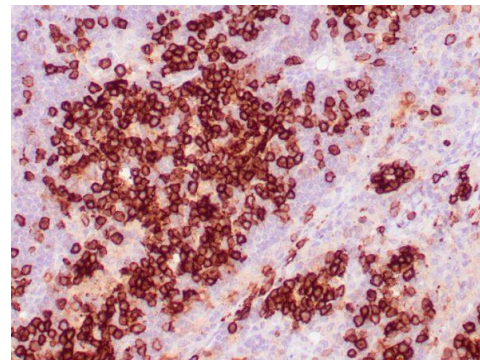

IHC analysis of formalin-fixed, paraffin-embedded human tonsil. Strong cytoplasmic and membranous staining using CD38/4328 at 2ug/ml in PBS for 30min RT. HIER: Tris/EDTA, pH9.0, 45min. 2 °: HRP-polymer, 30min. DAB, 5min.

Ordering Information and Current Pricing at www.scytek.com

Procedure:

- Tissue Section Pretreatment (Highly Recommended):** Staining of formalin fixed, paraffin embedded tissue sections is significantly enhanced by pretreatment with Tris-EDTA Solution (10x) pH 9.0 (ScyTek catalog# TES500) or Citrate Plus (10x) HIER Solution (ScyTek catalog# CPL500).
- Primary Antibody Incubation Time:** We suggest an incubation period of 30 minutes at room temperature. However, depending upon the fixation conditions and the staining system employed, optimal incubation should be determined by the user.
- Visualization:** For maximum staining intensity we recommend the "UltraTek HRP Anti-Polyvalent Lab Pack" (ScyTek catalog# UHP125, see IFU for instructions) combined with the "DAB Chromogen/Substrate Bulk Pack (High Contrast)" (ScyTek catalog# ACV500, see IFU for instructions).

Precautions:


- Contains Sodium Azide as a preservative (0.09% w/v).
- Do not pipette by mouth.
- Avoid contact of reagents and specimens with skin and mucous membranes.
- Avoid microbial contamination of reagents or increased nonspecific staining may occur.
- This product contains no hazardous material at a reportable concentration according to U.S. 29 CFR 1910.1200, OSHA Hazardous Communication Standard and EC Directive 91/155/EC.


References:

- Deaglio S, et. al. J Immunol. 1998;160(1):395-402.

Warranty:

No products or "Instructions For Use (IFU)" are to be construed as a recommendation for use in violation of any patents. We make no representations, warranties or assurances as to the accuracy or completeness of information provided on our IFU or website. Our warranty is limited to the actual price paid for the product. ScyTek Laboratories, Inc. is not liable for any property damage, personal injury, time or effort or economic loss caused by our products. Immunohistochemistry is a complex technique involving both histological and immunological detection methods. Tissue processing and handling prior to immunostaining can cause inconsistent results. Variations in fixation and embedding or the inherent nature of the tissue specimen may cause variations in results. Endogenous peroxidase activity or pseudoperoxidase activity in erythrocytes and endogenous biotin may cause non-specific staining depending on detection system used.

Storage: 2° C  8° C

 ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.